




NON ABSORBABLE SUTURE

SUTURES

# ARTHROV<sup>®</sup>LINK

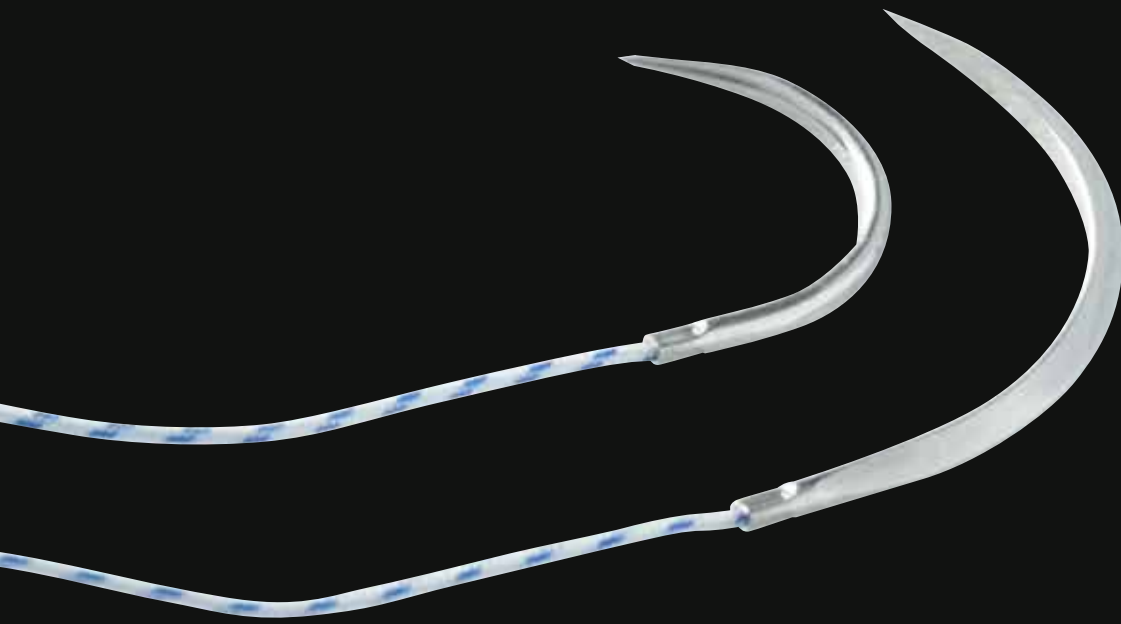

# NON ABSORBABLE SUTURE

## For tendon and ligament reinforcement and reinsertion.

**ARTHROLINK** is a non-absorbable braided sterile suture prepared from threads of ultra-high molecular weight, long-chain linear polyethylene (UHMWPE). It is indicated for repair or reinforcement of ligament, tendon and ligament reinsertion, tuberosities reinsertion, approximation and/or ligation of soft tissue.




### SUTURE NEEDLE Suture with double needle

**Type :**  
Suture USP2

**Size :**  
Length 90 cm

**Colors :**  
white/blue – white/black – white/green – blue.

**Needles :**  
Taper point 1/2 circle – length: 26 mm (transtendinous)  
Reverse cutting 3/8 circle – length: 40 mm (transosseous)

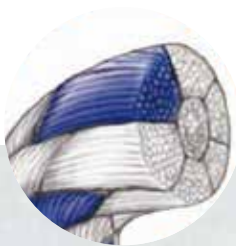

White / Blue



White / Black



White / Green



Blue

### FIBER STRUCTURE\*

Inter-locking core technology provides all fiber configurations with a sturdy core within center of the suture to act as a backbone for better knot tying and handling characteristics. This unique braid construction contributes to its superior strength. Its design is supple, soft for a better handling.





### SUTURE TAPE

#### Tape

**Type :**  
Suture USP2

**Size :**  
Tape length: 19 cm  
Tape width: 2 mm  
Suture Overall length: 100 cm

**Colors :**  
blue/blue – white/black

**Type :**  
Suture USP2

**Size :**  
Tape length: 60 cm  
Tape width: 1,6 mm  
Suture Overall length: 100 cm

**Colors :**  
white/green – white/blue



### SUTURE LOOP

#### Suture loop with needle

**Type :**  
Suture USP5

**Size :**  
Length 50 cm

**Needles & colors:**  
Triangular point 1/2 circle  
48 mm- blue

**Type :**  
Suture USP2


**Size :**  
Length 50 cm

**Needles & colors:**  
Taper point 1/2 circle  
26 mm- white/green  
Taper point 1/2 circle  
40 mm- white/blue

### TENSILE STRENGTH

**ARTHROLINK** offers excellent mechanical properties much stronger than polyester material and market leader. The braiding allows a better handling together with a high knot strength. Sutures remain stable after in-vivo implantation.

TENSILE STRENGTH\*

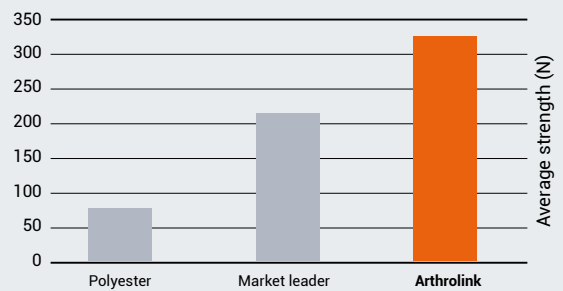

KNOT TENSILE STRENGTH\*



## SUTURES ADVANTAGES

- HIGH TENSILE STRENGTH
- SLIDING KNOT
- SOLIDITY AND SECURITY OF THE SUTURES
- BIOCOMPATIBLE
- 4 COLORS FOR A BETTER VISIBILITY



## STORAGE

**ARTHROLINK** must be stored in its original unopened packaging, in a clean and dry place, at ambient temperature.

## PRODUCTS

**ARTHROLINK** is packaged ready to use in a sealed pouch and sterilized by ethylene oxide. Single use. Do not re-sterilize.

**For any further information, please refer to the IFU.**

Designation		Colors	Reference
ARTHROLINK Suture Needle	USP2	white / blue	HSe-2C-VIMS
	USP2	white / black	HSe-2K-VIMS
	USP2	white / green	HSe-2G-VIMS
	USP2	blue	HS1008
ARTHROLINK Suture Loop	USP2	white / blue	HS955
	USP2	white / green	HS956
	USP5	Bleu	HS1042
ARTHROLINK Suture Tape	1,6mm	white / blue	HS1007
	1,6mm	white / green	HS1006
	2,0mm	blue / blue	HS957
	2,0mm	white / black	HS958



\*Manufacturer datas

Document design date : February 2020 - Range : Non absorbable suture -  
 Manufacturer : RIVERPOINT - Brand Name : ARTHROLINK - Medical Device  
 Class : III - CE Mark N° : 0459 User : Health professional - Indications :  
 Tendon and ligament reinforcement and reinsertion - Recommendations for  
 use: it is highly recommended to read the labels and the IFU.



VIMS SA  
 10 AVENUE DE FONTRÉAL - EUROCENTRE  
 31620 VILLENEUVE LES BOULOC - FRANCE  
 TÉL. +33 5 34 45 09 09 - FAX. +33 5 61 22 85 84

**VIMS-SYSTEM.COM**