




ALL SUTURE ANCHOR

SA  
SOFT ANCHOR

# ARTHRO **VIMS**<sup>®</sup> SOFT ANCHOR

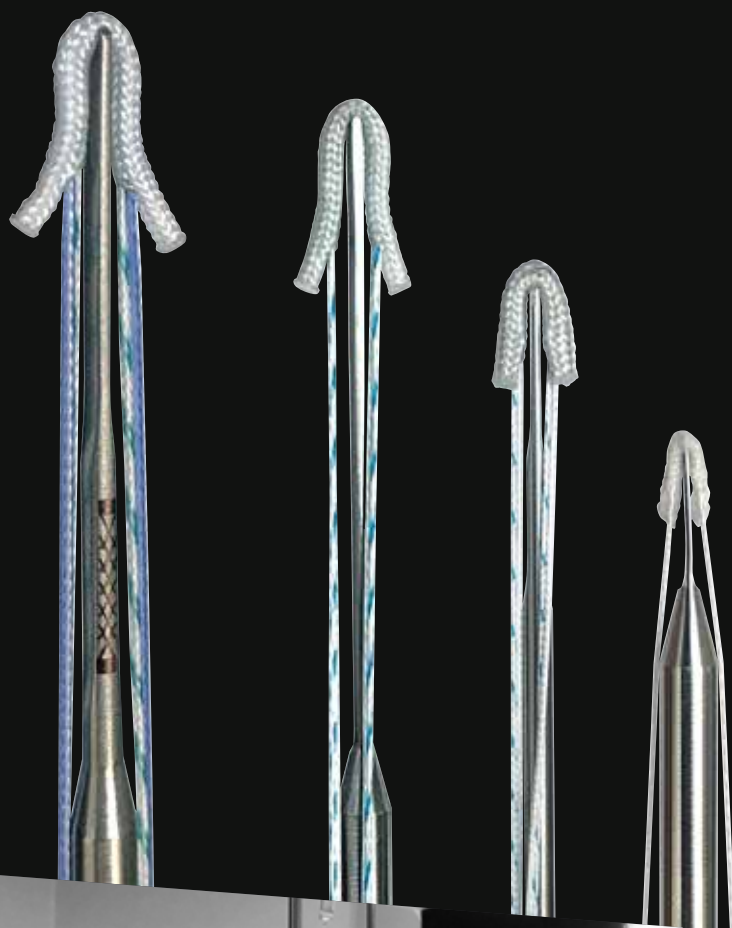


# ALL SUTURE ANCHOR

## Designed specially for the fixation of soft or prosthetic tissues to the bone

**ARTHROVIMS SOFT ANCHOR** are impacted anchor of small diameter exclusively composed of textile fibers.

The polyester body of the anchor is associated with one or two UHMWPE suture threads which are mounted on a disposable screwdriver allowing insertion into the bone.

The traction on the suture threads makes it possible to modify the shape of the body of the anchor which thus locks itself behind the cortical bone and ensures an effective anchorage.

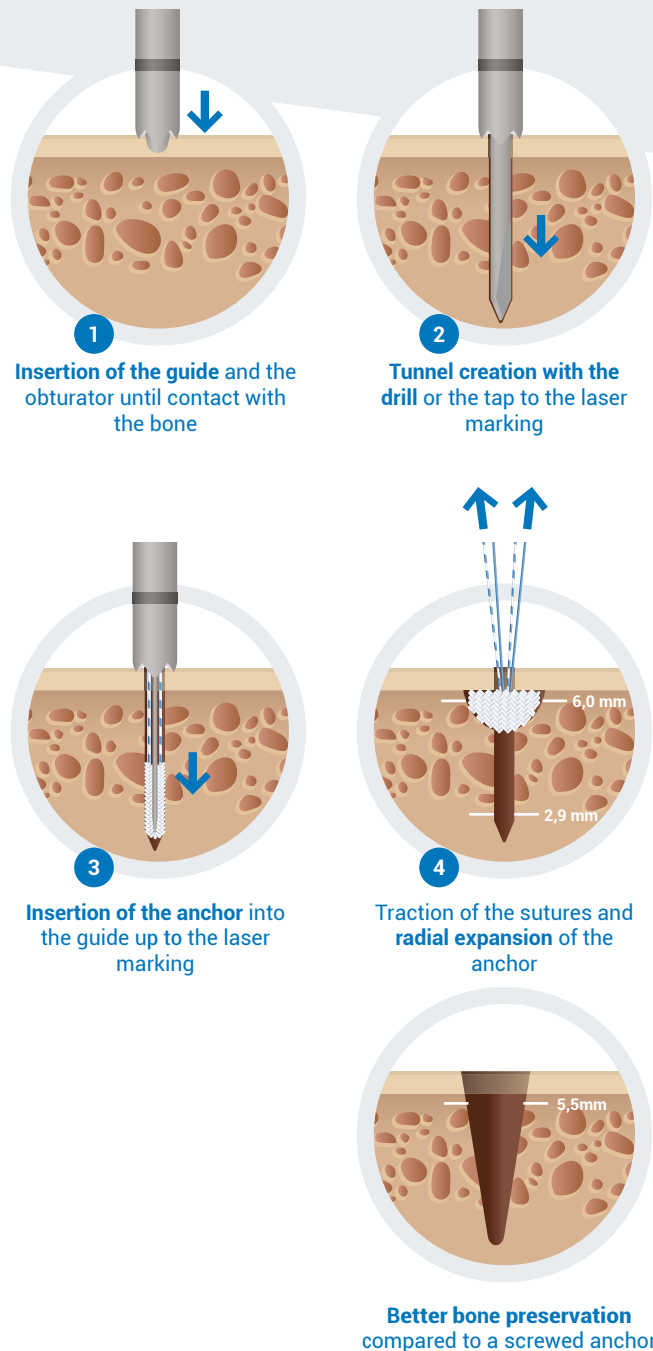

## ADVANTAGES

- Bone economy
- Easy to use
- No rigid material in case of revision or in the joint in case of tearing
- Radio-transparent
- Can be used under arthroscopy or mini open (with needles)
- Biocompatible
- High pull-out resistance
- Non-absorbable

**Beveled tip**  
for better impaction  
in bone tissue

**Laser marking**  
for insertion visualization





**SHOULDER**  
ROTATOR CUFF,  
BICEPS TENODESIS AND  
REINSERTION

**Ø2,9  
mm**



**SHOULDER**  
SLAP, BANKART  
INSTABILITY



**ANKLE**  
TENDON AND LIGAMENT  
REINSERTION

**Ø1,8  
mm**



**ANKLE**  
TENDON AND LIGAMENT  
REINSERTION



**SHOULDER**  
SLAP, BANKART  
INSTABILITY

**Ø1,5  
mm**



**HAND**  
TENDON AND LIGAMENT  
REINSERTION

**Ø1,2  
mm**



## SUTURE ADVANTAGES

ARTHROVIMS SOFT ANCHOR are associated with 2 or 3 USP2 sutures  
(made of very high molecular weight polyethylene UHMWPE)

HIGH TENSILE STRENGTH  
EASY REALIZATION OF KNOTS  
STRENGTH AND SAFETY OF SUTURES



## STORAGE

**ARTHROVIMS® SOFT ANCHOR** must be stored in its original unopened packaging, in a clean and dry place, at ambient temperature.

## PRODUCTS & INSTRUMENTATION

**ARTHROVIMS® SOFT ANCHOR** is packaged in a double pouch and sterilized by ethylene oxide. Single use.

**Do not re-sterilize.**

Instrumentation for arthroscopic procedures is available.

**ARTHROVIMS® SOFT ANCHOR** must be used with the corresponding instrumentation supplied by VIMS.

For any further information, please refer to the IFU.

DESIGNATION	SIZE	REFERENCE
ARTHROVIMS SOFT ANCHOR	2,9	VIMS-029
ARTHROVIMS SOFT ANCHOR with needles	2,9 needle	VIMS-029N
ARTHROVIMS SOFT ANCHOR	1,8	VIMS-018
ARTHROVIMS SOFT ANCHOR with needles	1,8 needle	VIMS-018N
ARTHROVIMS SOFT ANCHOR	1,5	VIMS-015
ARTHROVIMS SOFT ANCHOR with needles	1,5 needle	VIMS-015N
ARTHROVIMS SOFT ANCHOR with needles	1,2 needle	VIMS-012N

## INSTRUMENTATION

**ARTHROVIMS® SOFT ANCHOR** must be used with the specific instrumentation provided by VIMS.



CE 1984  
CLASS IIb / ANCHOR

CE  
CLASS I / INSTRUMENTATION

Document design date : February 2020 - Manufacturer : DORATEK - Range : All suture anchor - Brand Name : ARTHROVIMS SOFT ANCHOR, INSTRUMENTATION - User : Health professional - CE Mark N°: 1984 - Medical Device Class : IIb, I non stérile - Indications : Fixation of soft or prosthetic tissues to the bone - Recommendations for use: it is highly recommended to read the labels and the IFU.



VIMS SA  
10 AVENUE DE FONTRÉAL - EUROCENTRE  
31620 VILLENEUVE LES BOULOC - FRANCE  
TÉL. +33 5 34 45 09 09 - FAX. +33 5 61 22 85 84

VIMS-SYSTEM.COM