




ALL SUTURE ANCHOR

SOFT ANCHOR

SA

ARTHRO **VIMS**[®]

SOFT ANCHOR

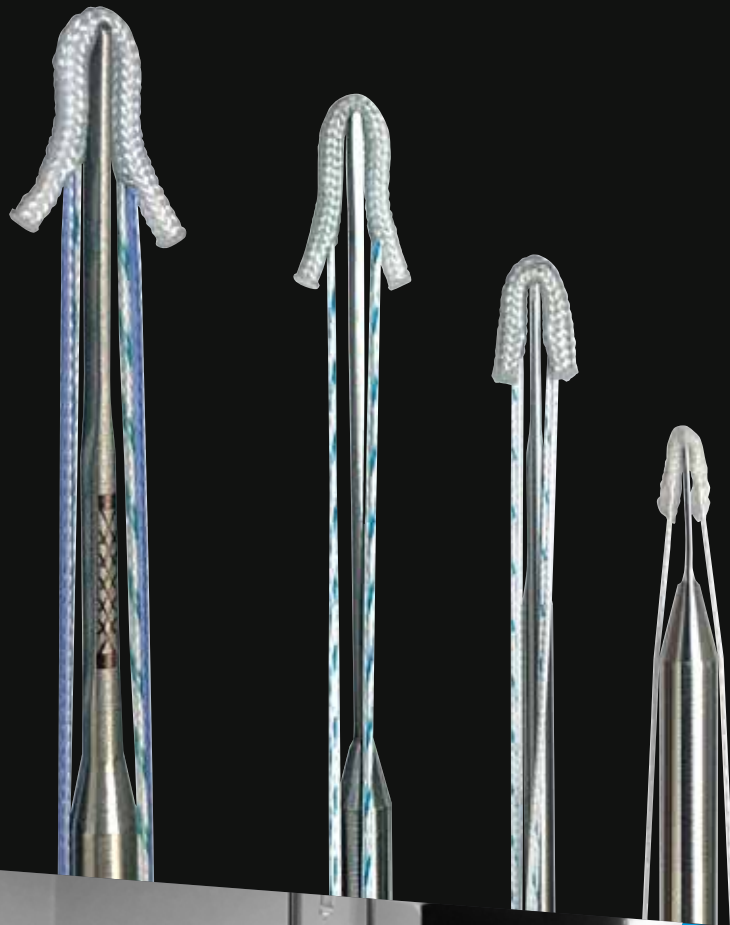




ALL SUTURE ANCHOR

It is specially designed for the fixation of soft or prosthetic tissues to the bone

ARTHROVIMS SOFT ANCHOR are impacted anchor of small diameter exclusively composed of textile fibers.

The polyester body of the anchor is associated with one or two UHMWPE suture threads which are mounted on a disposable screwdriver allowing insertion into the bone.

The traction on the suture threads makes it possible to modify the shape of the body of the anchor which thus locks itself behind the cortical bone and ensures an effective anchorage.

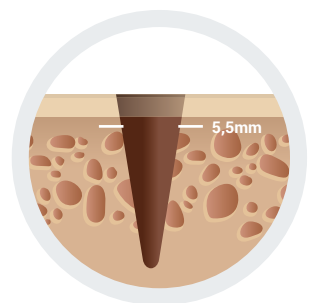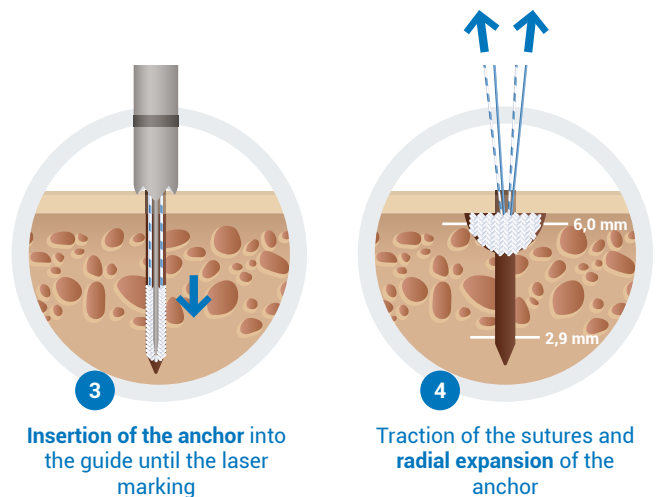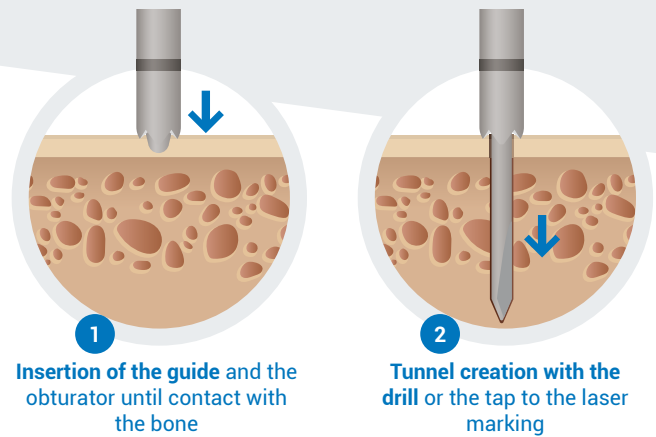

Bone economy compared to a screwed anchor

ADVANTAGES

- Bone economy
- Easy to use
- No rigid material in case of revision or in the joint in case of tearing
- Radio-transparent
- Can be used under arthroscopy or mini open (with needles)
- Biocompatible
- High pull-out resistance
- Non-absorbable

Beveled tip for better impaction in bone tissue

Laser marking for insertion visualization





SHOULDER
ROTATOR CUFF,
BICEPS TENODESIS AND
REINSERTION

Ø2,9
mm



SHOULDER
SLAP, BANKART
INSTABILITY



ANKLE
TENDON AND LIGAMENT
REINSERTION

Ø1,8
mm



ANKLE
TENDON AND LIGAMENT
REINSERTION



SHOULDER
SLAP, BANKART
INSTABILITY

Ø1,5
mm



HAND
TENDON AND LIGAMENT
REINSERTION

Ø1,2
mm



SUTURE ADVANTAGES

ARTHROVIMS SOFT ANCHOR are associated with 2 or 3 USP2 sutures
(made of very high molecular weight polyethylene UHMWPE)



HIGH TENSILE STRENGTH
EASY REALIZATION OF KNOTS
STRENGTH AND SAFETY OF SUTURES

STORAGE

ARTHROVIMS® SOFT ANCHOR must be stored in its original unopened packaging, in a clean and dry place, at ambient temperature.

PRODUCTS & INSTRUMENTATION

ARTHROVIMS® SOFT ANCHOR is packaged in a double pouch and sterilized by ethylene oxide. Instrumentation for arthroscopic procedures is available. Single use. **Do not re-sterilize.**

ARTHROVIMS® SOFT ANCHOR must be used with the corresponding instrumentation supplied by VIMS. For any further information, please refer to the IFU.

DESIGNATION	SIZE	REFERENCE
ARTHROVIMS SOFT ANCHOR (2 sutures)	2,9	VIMS-029
ARTHROVIMS SOFT ANCHOR (2 sutures with needles)	2,9	VIMS-029N
ARTHROVIMS SOFT ANCHOR (2 sutures)	1,8	VIMS-018
ARTHROVIMS SOFT ANCHOR (2 sutures with needles)	1,8	VIMS-018N
ARTHROVIMS SOFT ANCHOR (1 suture)	1,5	VIMS-015
ARTHROVIMS SOFT ANCHOR (1 suture with needle)	1,5	VIMS-015N
ARTHROVIMS SOFT ANCHOR (1 suture with needle)	1,2	VIMS-012N
ARTHROVIMS SOFT ANCHOR (1 suture + 1 tape blue)	2,9	VIMS-29FTB
ARTHROVIMS SOFT ANCHOR (1 suture + 1 tape white/blue)	2,9	VIMS-029FTWB

INSTRUMENTATION

ARTHROVIMS® SOFT ANCHOR must be used with the specific instrumentation provided by VIMS.



CE 1984
CLASSE IIb / ANCHOR

CE
CLASSE I / INSTRUMENTATION

Document design date : Mai 2019 - Manufacturer : DORATEK - Range : All suture anchor - Brand Name : ARTHROVIMS SOFT ANCHOR, INSTRUMENTATION - User : Health professional - CE Mark N°: 1984 - Medical Device Class : IIb, I non stérile - Indications : Fixation of soft or prosthetic tissues to the bone - Recommendations for use: it is highly recommended to read the labels and the IFU.

VIMS

VIMS SA
10 AVENUE DE FONTRÉAL - EUROCENTRE
31620 VILLENEUVE LES BOULOC - FRANCE
TÉL. +33 5 34 45 09 09 - FAX. +33 5 61 22 85 84

VIMS-SYSTEM.COM