




BUTTON TENDINOUS OR LIGAMENOUS
FIXATION SYSTEM

BUTTON

ARTHRO **VIMS**[®]

BUTTON

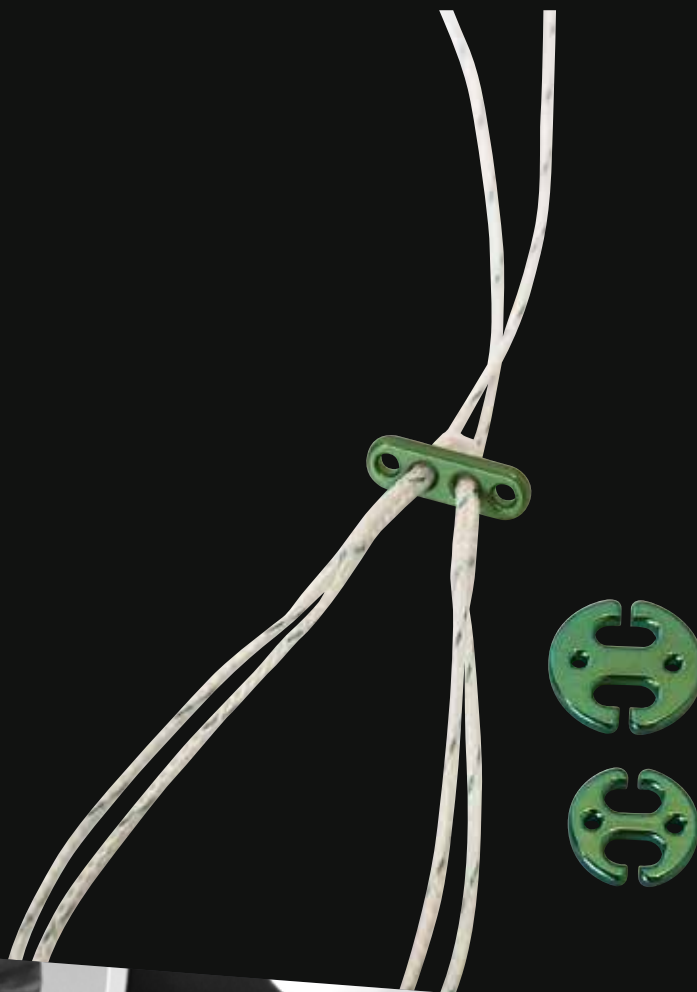

V | **SPORT
MEDICINE**
ADVANCED TECHNOLOGY

ARTHROVIMS BUTTON

Tendinous or ligamentous fixation system

VIMS offers a wide range of fixation devices dedicated to minimal invasive anterior cruciate ligament reconstruction surgeries (button with pre-measured loop size or adjustable version), repair of acromioclavicular joint dislocation, fracture of the lateral quarter of the clavicle and Latarjet procedures.



SHOULDER



KNEE



BUTTON LOOP

- Pre-measured loop size button is made of :
- 1 titanium button (dimensions 12x4 mm)
 - Polyester PET braided pre-measured loop
 - 2 pre mounted UHMWPE sutures (green and white)
 - 1 extension titanium plate (optional)

Indications

Tendon or ligament fixation in reconstruction surgeries especially for ACL repair.

Loop sizes

10 sizes are available and meet requirements of conventional ligamentoplasty surgeries (15, 20, 25, 30, 35, 40, 45, 50, 55, 60mm)





ADJUSTABLE BUTTON LOOP (KNEE)

Adjustable button is made of :

- 1 titanium button (dimensions 12x4 mm)
- 1 UHMWPE braided double loop
- 2 pre mounted UHMWPE pulling wire (white and white/black)
- 1 extension titanium plate (optional)

Indications

Tendon or ligament fixation in reconstruction surgeries and especially for ACL ligamentoplasty.



ADJUSTABLE BUTTON LOOP (SHOULDER)

System is made of two devices :

- 1 adjustable button loop external (titanium plate, UHMWPE double loop and pulling wires)
- 1 or 2 open titanium buttons in 8 mm or 10 mm size

Indications

Repair of acromium clavicular joint dislocation, lateral quarter of the clavicle fracture, Latarjet procedures.



STORAGE

ARTHROVIMS® BUTTON devices must be stored in original sealed packaging, in dry and clean environment, at ambient temperature.

For any further information, please read the instructions for use.

PRODUCTS

ARTHROVIMS® BUTTON devices are delivered ready to use, sterile and packed in double sealed pouch.

Sterilization method is EO according to the Norm ISO 11135

Single use. **Do not re-use or re-process.**

Instruments for ACL/PCL ligamentoplasty, acromio-clavicular dislocation, Latarjet are **available upon request**. Do not hesitate to contact us for more details.

DESIGNATION	SIZE/ COLOR	REFERENCE
ADJUSTABLE BUTTON LOOP external	White/Blue	VSDL-2001
ADJUSTABLE BUTTON LOOP external	White/Green	VSDL-2002
ADJUSTABLE BUTTON LOOP internal	White/Blue	VLFT-1000
TITANIUM SHOULDER BUTTON	8 mm Blue	VDL-2000B
TITANIUM SHOULDER BUTTON	8 mm Gren	VDL-2000G
TITANIUM SHOULDER BUTTON	10 mm Blue	VDL-2010B
TITANIUM SHOULDER BUTTON	10 mm Green	VDL-2010G
BUTTON LOOP	15 mm	VBL-15
	20 mm	VBL-20
	25 mm	VBL-25
	30 mm	VBL-30
	35 mm	VBL-35
	40 mm	VBL-40
	45 mm	VBL-45
	50 mm	VBL-50
	55 mm	VBL-55
	60 mm	VBL-60
BUTTON LOOP EXTENSION PLATE	Blue	VEP-01
LATARJET STANDARD KIT ¹	8 mm	VLSK-06
LATARJET SMALL SHOULDER KIT ²	10 mm	VLSS-03

¹ **LATARJET STANDARD KIT** : 2 adjustable fixation systems, 2 buttons size 8 mm, 2 sutures UHMWPE.

² **LATARJET SMALL SHOULDER KIT** : 1 adjustable fixation system, 1 button size 10 mm, 1 suture UHMWPE.

