




ALL SUTURE ANCHOR

SA /
SOFT ANCHOR

ARTHRO **VIMS**®

SOFT ANCHOR

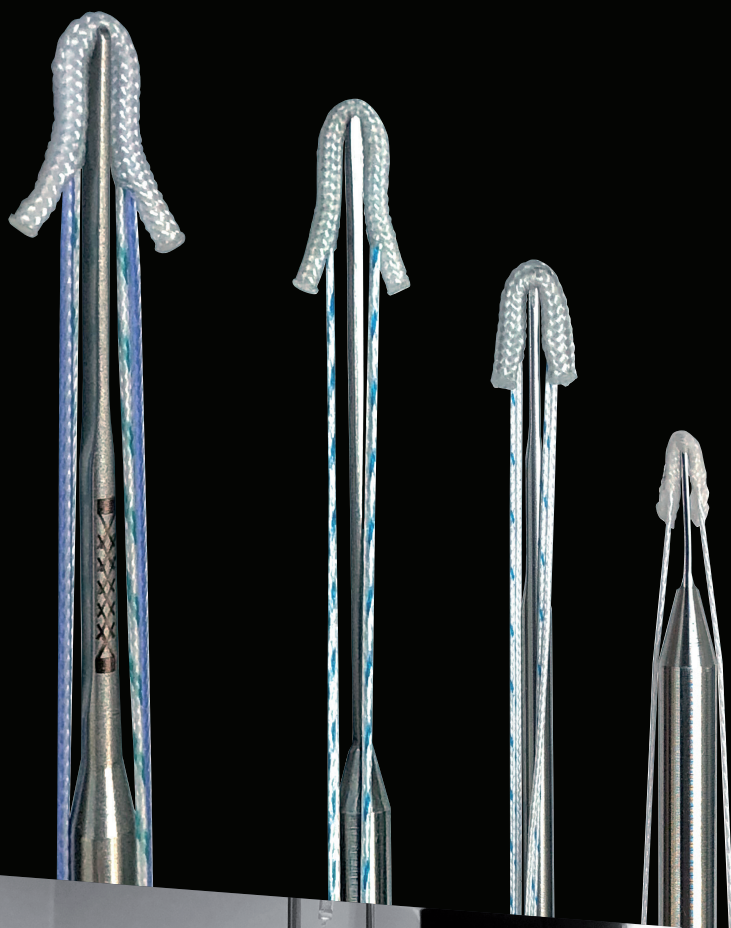


ALL SUTURE ANCHOR

It is specially designed for the fixation of soft or prosthetic tissues to the bone

ARTHROVIMS SOFT ANCHOR are impacted anchor of small diameter exclusively composed of textile fibers.

The polyester body of the anchor is associated with one or two UHMWPE suture threads which are mounted on a disposable screwdriver allowing insertion into the bone.

The traction on the suture threads makes it possible to modify the shape of the body of the anchor which thus locks itself behind the cortical bone and ensures an effective anchorage.

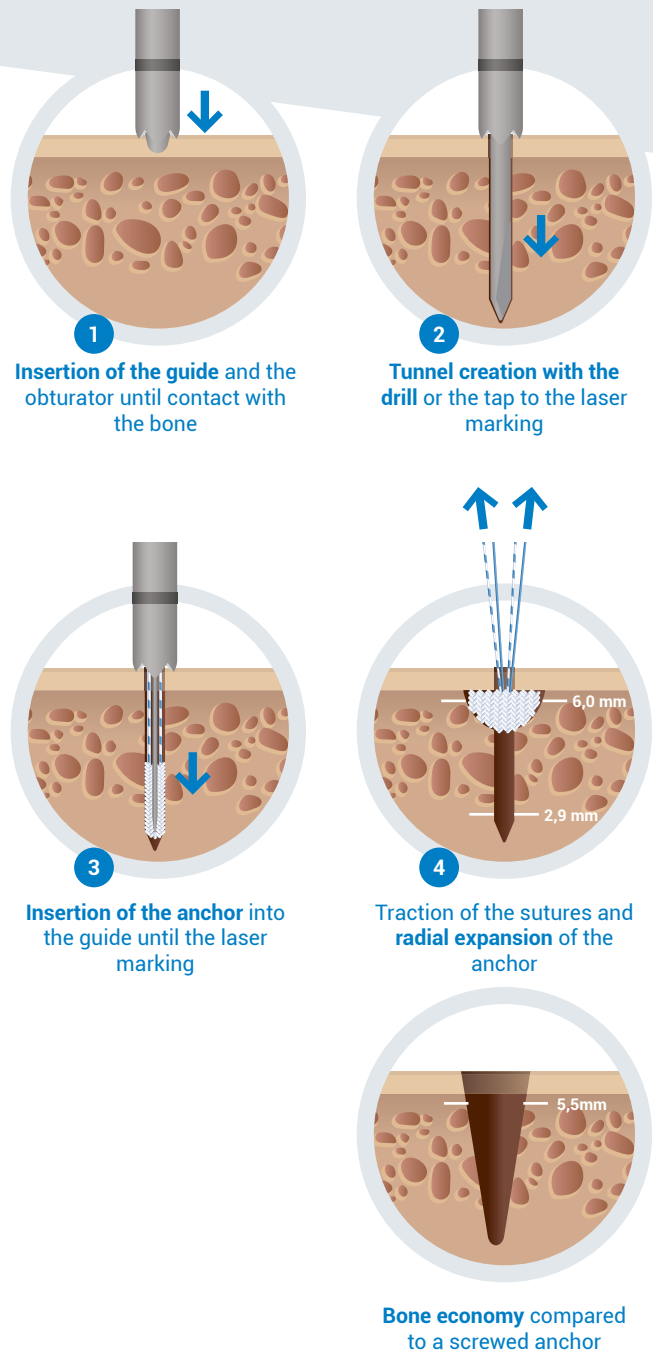

ADVANTAGES

- Bone economy
- Easy to use
- No rigid material in case of revision or in the joint in case of tearing
- Radio-transparent
- Can be used under arthroscopy or mini open (with needles)
- Biocompatible
- High pull-out resistance
- Non-absorbable

Beveled tip
for better impaction
in bone tissue


Laser marking
for insertion visualization





SHOULDER
ROTATOR CUFF,
BICEPS TENODESIS AND
REINSERTION

**Ø2,9
mm**

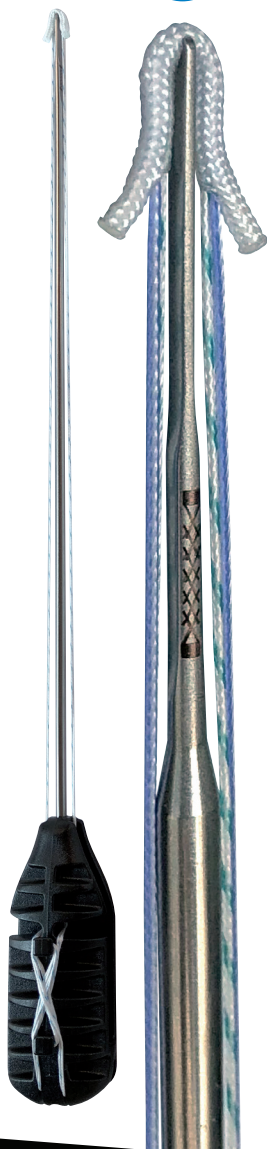


SHOULDER
SLAP, BANKART
INSTABILITY



ANKLE
TENDON AND LIGAMENT
REINSERTION

**Ø1,8
mm**

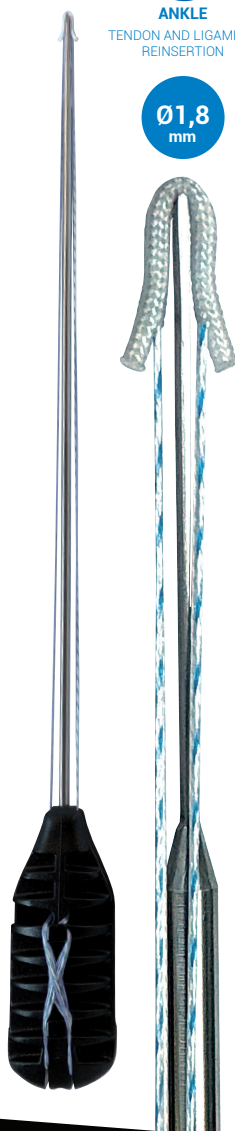

ANKLE
TENDON AND LIGAMENT
REINSERTION

**Ø1,5
mm**



HAND
TENDON AND LIGAMENT
REINSERTION


**Ø1,2
mm**



SUTURE ADVANTAGES

ARTHROVIMS SOFT ANCHOR are associated with 2 or 3 USP2 sutures
(made of very high molecular weight polyethylene UHMWPE)

HIGH TENSILE STRENGTH
EASY REALIZATION OF KNOTS
STRENGTH AND SAFETY OF SUTURES

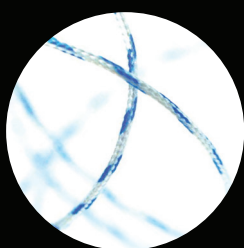

STORAGE

ARTHROVIMS® SOFT ANCHOR must be stored in its original unopened packaging, in a clean and dry place, at ambient temperature.

PRODUCTS & INSTRUMENTATION


ARTHROVIMS® SOFT ANCHOR is packaged in a double pouch and sterilized by ethylene oxide. Instrumentation for arthroscopic procedures is available. Single use. **Do not re-sterilize.**

ARTHROVIMS® SOFT ANCHOR must be used with the corresponding instrumentation supplied by VIMS. For any further information, please refer to the IFU.

DESIGNATION	SIZE	REFERENCE
ARTHROVIMS SOFT ANCHOR (2 sutures)	2,9	VIMS-029
ARTHROVIMS SOFT ANCHOR (2 sutures with needles)	2,9	VIMS-029N
ARTHROVIMS SOFT ANCHOR (2 sutures)	1,8	VIMS-018
ARTHROVIMS SOFT ANCHOR (2 sutures with needles)	1,8	VIMS-018N
ARTHROVIMS SOFT ANCHOR (1 suture with needle)	1,5	VIMS-015N
ARTHROVIMS SOFT ANCHOR (1 suture with needle)	1,2	VIMS-012N
ARTHROVIMS SOFT ANCHOR (1 suture + 1 tape blue)	2,9	VIMS-29FTB
ARTHROVIMS SOFT ANCHOR (1 suture + 1 tape white/blue)	2,9	VIMS-029FTWB

INSTRUMENTATION

ARTHROVIMS® SOFT ANCHOR must be used with the specific instrumentation provided by VIMS.

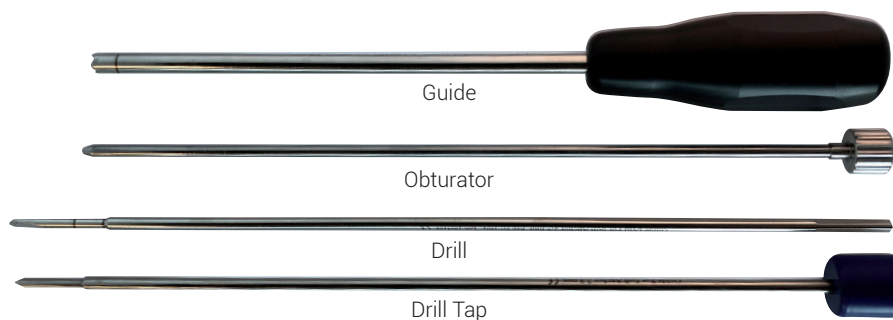

CE 2292
CLASSE IIb / ANCHOR

CE
CLASSE I / INSTRUMENTATION