


ARTHROV[®]LINK[™]

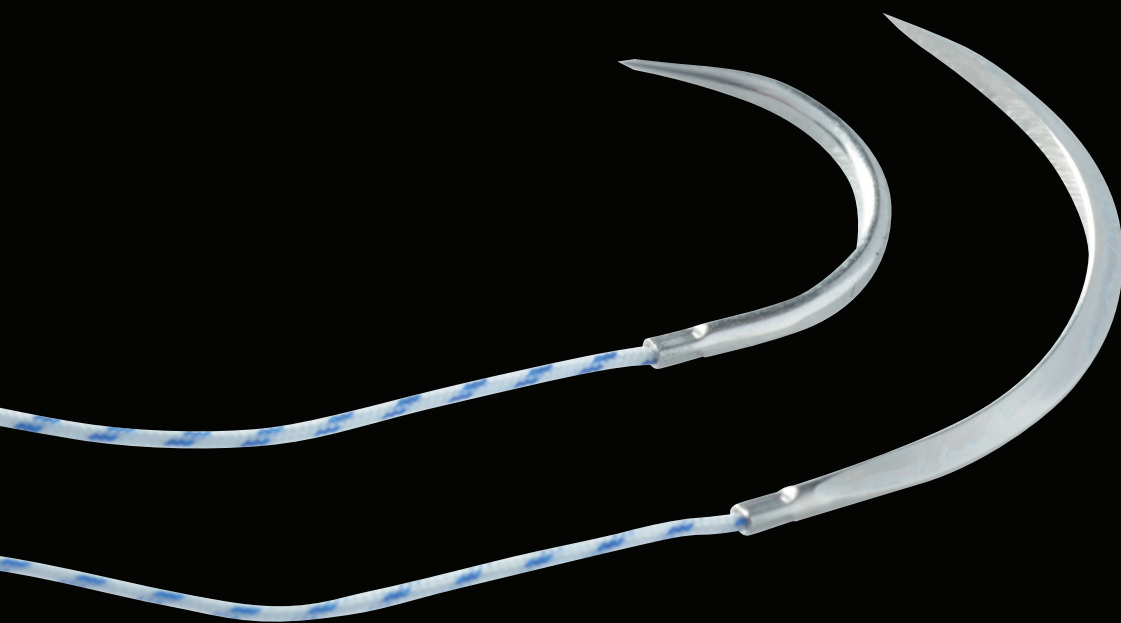

NON ABSORBABLE SUTURE

For tendon and ligament reinforcement and reinsertion

ARTHROLINK is a non-absorbable braided sterile suture prepared from threads of ultra-high molecular weight, long-chain linear polyethylene (UHMWPE). It is indicated for repair or reinforcement of ligament, tendon and ligament reinsertion, tuberosities reinsertion, approximation and/or ligation of soft tissue.



SHOULDER



HIP




KNEE



ANKLE



White / Blue

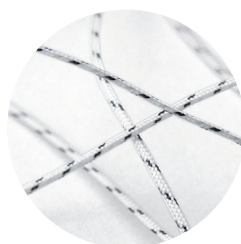

White / Black



White / Green



Blue

FIBER STRUCTURE*

Inter-locking core technology provides all fiber configurations with a sturdy core within center of the suture to act as a backbone for better knot tying and handling characteristics. This unique braid construction contributes to its superior strength. Its design is supple, soft for a better handling.



SUTURE NEEDLE Suture with double needle

Type :
Suture USP2
Suture USP5

Size :
Length 90 cm

Colors :
white/blue – white/black – white/green – blue.

Needles :
Taper point 1/2 circle – length: 26 mm (transtendinous)
Reverse cutting 3/8 circle – length: 40 mm (transosseous)





SUTURE TAPE

Tape



Type	Tape length	Tape width	Suture length	Colors	Needles
USP2	19 cm	2,0 mm	100 cm	Blue/Blue or White/Black	No
USP2	21 cm	2,0 mm	100 cm	White/Blue	No
USP2	60 cm	1,6 mm	100 cm	White/Green or White/Blue	No
USP2	80 cm	1,6 mm	100 cm	Blue/Black or Black	Taper point 1/2 circle 26 mm (transtendinous) Reverse cutting 3/8 circle - 40mm (transosseous)

SUTURE LOOP

Name	Type	Length	Colors	Needles
Suture loop with needles	USP0	50cm	Black	Straight round needle 60mm
Suture loop with needles	USP5	50 cm	Blue or White/Black	Triangular point 1/2 circle 48 mm
Suture loop with needles	USP2	50 cm	White/Green	Taper point 1/2 circle 26 mm
Suture loop with needles	USP2	50 cm	White/Blue	Taper point 1/2 circle 48 mm
Suture Shuttle Loop	USP2	66 cm	Blue	No

TENSILE STRENGTH

ARTHROLINK offers excellent mechanical properties much stronger than polyester material and market leader. The braiding allows a better handling together with a high knot strength. Sutures remain stable after in-vivo implantation.

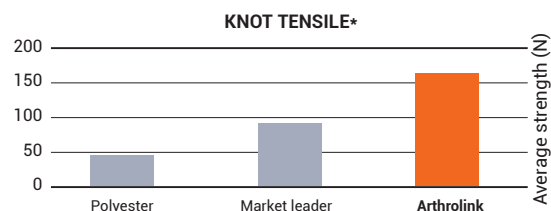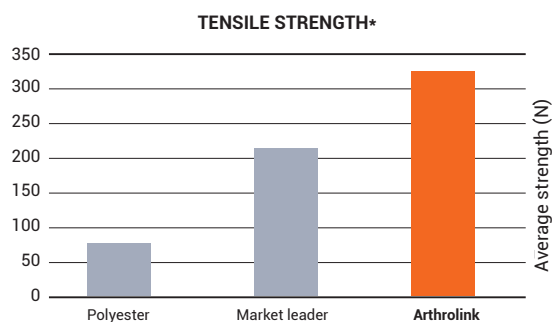

SUTURES ADVANTAGES

- HIGH TENSILE STRENGTH
- KNOT SLIDING
- SOLIDITY AND SECURITY OF THE SUTURES
- BIOCOMPATIBLE
- DIFFERENT COLORS FOR A BETTER VISIBILITY



STORAGE

ARTHROLINK sutures must be stored in their original, unopened packaging, in a clean and dry place at ambient temperature.

PRODUCTS

ARTHROLINK sutures are packaged in a ready-to-use pouch and sterilized by ethylene oxide.
Single use only. Do not re-sterilize.
For any further information, please read the IFU and/or contact us.

Designation	Size	Colors	Reference
ARTHROLINK Suture Needle	USP2	white / blue	HSe-2C-VIMS
	USP2	white / black	HSe-2K-VIMS
	USP2	white / green	HSe-2G-VIMS
	USP2	blue	HS1008
	USP5	white / blue	HS1043
	USP5	white / black	HS1044
	USP5	white / green	HS1045
	USP5	blue	HS1046
ARTHROLINK Suture Loop	USP2	white / blue needled	HS955
	USP2	white / green needled	HS956
	USP5	blue needled	HS1042
	USP5	white / black needled	HS1041
	USP0	black / straight needle	HS345
ARTHROLINK Suture Shuttle Loop	USP2	blue	HS1063
ARTHROLINK Suture Tape	1,6mm	white / blue	HS1007
	1,6mm	white / green	HS1006
	1,6mm	blue / black needled	HS1061
	1,6mm	black needled	HS1062
	2,0mm	blue / blue	HS957
	2,0mm	white / black	HS958
	2,0mm	white / blue	HS216



*Manufacturer datas

Document design date : March 2021 - Range : Non absorbable suture -
Manufacturer : RIVERPOINT - Brand Name : ARTHROLINK - Medical Device
Class : IIb - CE Mark N° : 0459 User : Health professional - Indications :
Tendon and ligament reinforcement and reinsertion - Recommendations for
use: it is highly recommended to read the labels and the IFU.
BR_ARTHROLINK_EN_V02



VIMS SA
10 AVENUE DE FONTRÉAL - EUROCENTRE
31620 VILLENEUVE LES BOULOC - FRANCE
TÉL. +33 5 34 45 09 09 - FAX. +33 5 61 22 85 84

VIMS-SYSTEM.COM

A  ZIMMER BIOMET company