

ARTHRO**VIMS**[®]
ANCHOR PK



PEEK SUTURE ANCHOR

**It is specially designed
for the fixation of soft or
prosthetic tissues to the bone.**

ARTHROVIMS® ANCHOR PEEK

are screwed or impacted anchors associated with two or three UHMWPE suture threads (USP2) and/or tapes 1,4mm and mounted on a disposable screwdriver for insertion and locking into the bone.

These anchors are composed of a thermoplastic polymer : the PEEK (polyetheretherketone) VESTAKEEP considered as one of the most successful in the world with more than 15 years of clinical experience in orthopedics.



SHOULDER
ROTATOR CUFF, BICEPS
TENODESIS AND
REINSERTION

DIAMETER
(mm)

Ø4,5
mm

Ø5,5
mm

Ø6,5
mm



ARTHROVIMS PK ANCHOR SCREWED TO CORTICAL LOCK

ARTHROVIMS PK anchor is a fully threaded suture anchor featuring dual threads to maximize cortical and cancellous fixation, has a rounded tip to protect the sutures and to facilitate the insertion.

Particularly suitable for repairing rotator cuff and associated pathologies, this anchor is available with two or three sutures and/or tapes.

Double cancellous and cortical thread

- Maximize fixation strength by increasing the interaction between cortical and cancellous bone
- Greater torsional torque when inserted
- Maximum tensile strength, therefore pull-out

Optimal internal foot print of the screwdriver

- More important contact surface between the anchor and the screwdriver for better torque and guided insertion


4,5 - 5,5 - 6,5 mm

2 or 3 sutures - 1 suture/1 tape - 2 tapes

- Adaptation to multiple surgical techniques

PEEK ADVANTAGES

- High mechanical properties
- Excellent biocompatibility
- Chemically inert material
- Radio-transparent
- Elastic modulus close to bone tissue : more homogeneous constraints transfer and reduction of stress zones

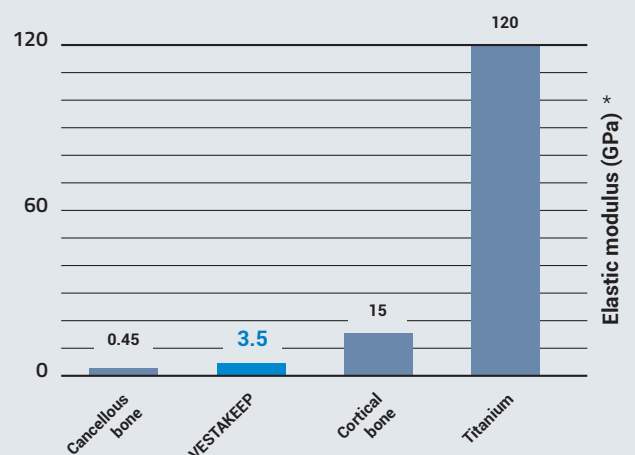

ARTHROVIMS LRP & LRT

IMPACTED SCREWED ANCHOR WITHOUT KNOT

ARTHROVIMS LRP et **LRT** are fully threaded, knotless anchor designed to be used with sutures or tapes for rotator cuff repair using the "bridge" techniques.

The "knotless" technique consists of passing sutures or tapes through tissue, loading them through the anchor eyelet which is then inserted into the bone socket.

This technique eliminates possible complications caused by knots compared to other conventional anchors.



SHOULDER
ROTATOR CUFF



ANKLE
ACHILLE TENDON

DIAMETER
(mm)

Ø4,75
mm

Ø5,5
mm



ARTHROVIMS MICRO-PK MINI-PK

IMPACT ANCHOR

ARTHROVIMS MICRO-PK and **MINI-PK** suture anchors are available in many configurations: several sizes, with one or two sutures, with or without needle or knotless. They present an eyelet which provides superior abrasion resistance due to PEEK's low coefficient of friction.

The use of **MICRO-PK** and **MINI-PK** anchors implements a simple and reproducible technique using a guide for the drilling and for the impactation of the anchor



SHOULDER
SLAP, BANKART
INSTABILITY



ANKLE
TENDON AND LIGAMENT
REINSERTION



HAND
REINSERTION

DIAMETER
(mm)

Ø1,8
mm

Ø2,0
mm

Ø2,8
mm

Ø3,0
mm

Ø3,5
mm



SHOULDER
BICEPS TENODESIS



ANKLE
ACHILLE TENDON

DIAMETER
(mm)

Ø5,0
mm

Ø6,0
mm

Ø7,0
mm



ARTHROVIMS BT

SCREWED ANCHOR WITHOUT KNOTS

ARTHROVIMS BT tenodesis implant are designed for all-arthroscopic or mini-open proximal biceps tenodesis. This system reduces the operating time by its simplicity and efficiency.


The implant feature a unique PEEK eyelet which is used to steer the biceps tendon to the bottom of the bone socket without need to externalize or whipstitch the tendon. It is not necessary to predetermine tendon length, as position and tension are easily set and visualized prior to final anchor insertion. Strong fixation is obtained by simply advancing the pre-loaded **ARTHROVIMS BT** anchor.

Drills and a graft calibration plate are available.

SUTURES ADVANTAGES

ARTHROVIMS ANCHOR PEEK are associated with 2 or 3 USP2 sutures and/or tapes of 1,4mm (made of very high molecular weight polyethylene UHMWPE)

HIGH TENSILE STRENGTH
EASY REALIZATION OF KNOTS
STRENGTH AND SAFETY OF SUTURES

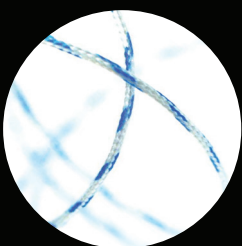

STORAGE

ARTHROVIMS® ANCHOR PEEK must be stored in its original unopened packaging, in a clean and dry place, at ambient temperature.

PRODUCTS & INSTRUMENTATION

ARTHROVIMS® ANCHOR PEEK is packaged in double pouch, and sterilized by ethylene oxide. Single use. Do not re-sterilize.

Instrumentation for for arthroscopic and mini-open surgery is available **ARTHROVIMS® ANCHOR PEEK** must be used with the corresponding instruments supplied by VIMS. For any further information, please refer to the IFU and/or contact us.

ANCHORS	SIZE	REFERENCES
ROTATOR CUFF		
ARTHROVIMS PK - 2 sutures	4,5	VIMS-9141.45F
ARTHROVIMS PK - 3 sutures	4,5	VIMS-9145.45F
ARTHROVIMS PK - 2 sutures	5,5	VIMS-9090.55F
ARTHROVIMS PK - 3 sutures	5,5	VIMS-9091.55F
ARTHROVIMS PK - 1 suture blue + 1 tape white/blue	4,5	VIMS-9180.45FT
ARTHROVIMS PK - 1 suture white/blue + 1 tape blue	4,5	VIMS-9181.45FT
ARTHROVIMS PK - 1 suture bleu + 1 tape black	5,5	VIMS-9182.55FT
ARTHROVIMS PK - 1 suture white/blue + 1 tape blue	5,5	VIMS-9183.55FT
ARTHROVIMS PK - 2 tapes white/blue & black + needles	4,5	VIMS-9222.45T
ARTHROVIMS PK - 2 tapes white/black & blue + needles	4,5	VIMS-9223.45T
ARTHROVIMS PK - 2 tapes white/blue & black + needles	5,5	VIMS-9224.55T
ARTHROVIMS PK - 2 tapes white/black & blue + needles	5,5	VIMS-9225.55T
ARTHROVIMS LRP - peek tip	5,5	VIMS-9092.55F
ARTHROVIMS LRP - peek tip	4,75	VIMS-9142.475F
ARTHROVIMS LRT - titanium tip	5,5	VIMS-9136.55F
ARTHROVIMS LRT - titanium tip	4,75	VIMS-9143.475F
SLAP, BANKART, INSTABILITY		
ARTHROVIMS MICRO PK - 1 suture + needles	1,8	VIMS-9139.18F
ARTHROVIMS MICRO PK - 1 suture + needles	2	VIMS-9140.20F
ARTHROVIMS MINI PK - 1 suture	2,8	VIMS-9086.28F
ARTHROVIMS MINI PK - 2 sutures	3	VIMS-9088.30F
ARTHROVIMS MINI PK - 2 sutures	3,5	VIMS-9089.35F
ARTHROVIMS MINI PK - Knotless	2,8	VIMS-9087.28F
ARTHROVIMS MINI PK - 2 sutures + needles	3	VIMS-9186.30F
BICEPS TENODESIS		
ARTHROVIMS BT	5	VIMS-9093.50F
ARTHROVIMS BT	6	VIMS-9094.60F
ARTHROVIMS BT	7	VIMS-9095.70F



Document design date : March 2021 - Manufacturer : Chetan Meditech - Range : PEEK suture anchor - Brand Name : ARTHROVIMS Anchor PEEK, Anchor PK, Anchor LRP, Anchor LRT, Anchor Micro PK, Anchor Mini PK, Anchor BT - User : Health professional - CE Mark N°: 2265 - Medical Device Class : IIb - Indications : fixation of soft or prosthetic tissues to the bone - Recommendations for use: it is highly recommended to read the labels and the IFU.
BR_ARTHROVIMS_ANCHOR_PEEK_EN_V02

*Datas provided by the manufacturer



VIMS SA
10 AVENUE DE FONTRÉAL - EUROCENTRE
31620 VILLENEUVE LES BOULOC - FRANCE
TÉL. +33 5 34 45 09 09 - FAX. +33 5 61 22 85 84

VIMS-SYSTEM.COM

A  ZIMMER BIOMET company