


ARTHRO**VIMS**[®]
SOFT ANCHOR

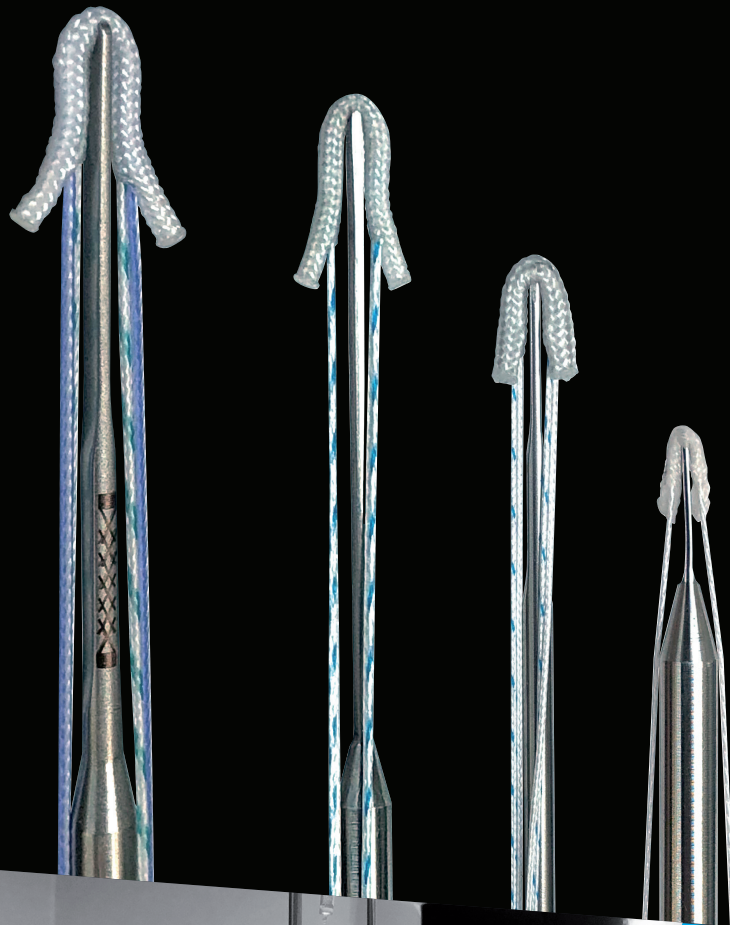




ALL SUTURE ANCHOR

It is specially designed for the fixation of soft or prosthetic tissues to the bone

ARTHROVIMS SOFT ANCHOR are impacted anchor of small diameter exclusively composed of textile fibers.

The polyester body of the anchor is associated with one or two UHMWPE suture threads, and/or tape, which are mounted on a disposable inserter allowing insertion into the bone.

The traction on the suture threads makes it possible to modify the shape of the body of the anchor which thus locks itself behind the cortical bone and ensures an effective anchorage.

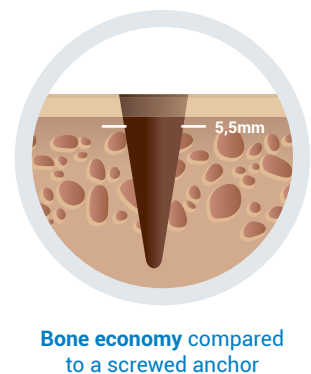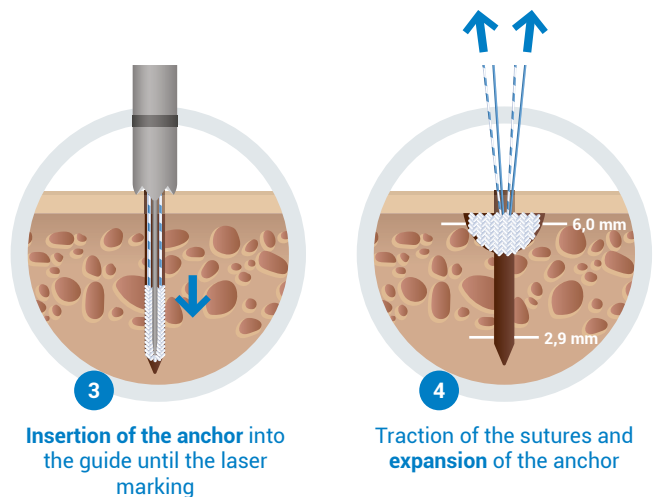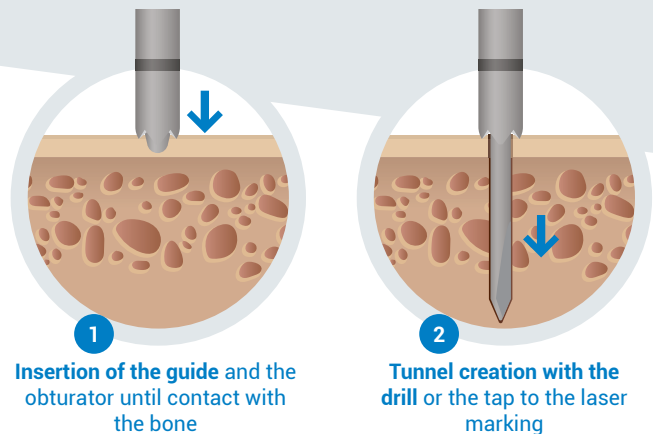

ADVANTAGES

- Bone economy
- Easy to use
- Radio-transparent
- Can be used under arthroscopy or mini open (with needles)
- Biocompatible
- High pull-out resistance
- Non-absorbable

Beveled tip
for better impaction
in bone tissue


Laser marking
for insertion visualization





SHOULDER
ROTATOR CUFF,
BICEPS TENODESIS AND
REINSERTION

Ø2,9
mm

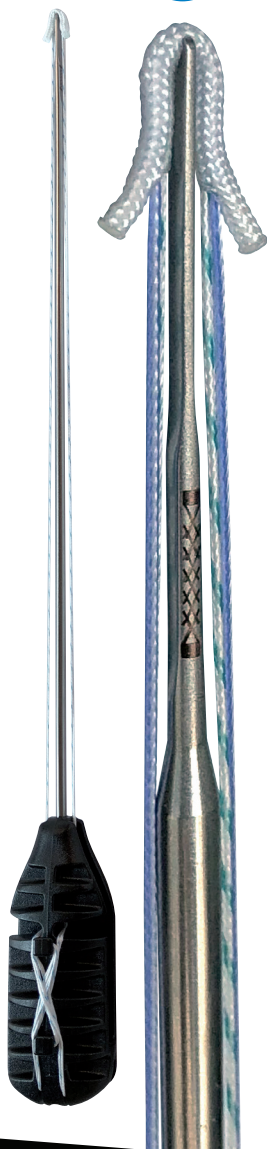


SHOULDER
SLAP, BANKART
INSTABILITY



ANKLE
TENDON AND LIGAMENT
REINSERTION

Ø1,8
mm

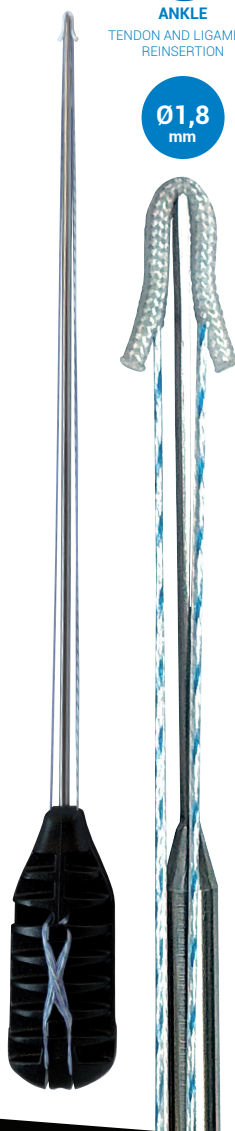

ANKLE
TENDON AND LIGAMENT
REINSERTION

Ø1,5
mm



HAND
TENDON AND LIGAMENT
REINSERTION


Ø1,2
mm



SUTURE ADVANTAGES

ARTHROVIMS SOFT ANCHOR are associated with 2 or 3 USP2 sutures
or tapes (made of very high molecular weight polyethylene UHMWPE)

HIGH TENSILE STRENGTH
EASY REALIZATION OF KNOTS
STRENGTH AND SAFETY OF SUTURES

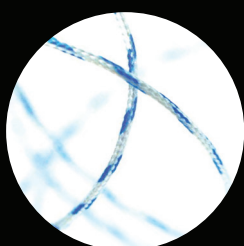

STORAGE

ARTHROVIMS® SOFT ANCHOR must be stored in its original unopened packaging, in a clean and dry place, at ambient temperature.


PRODUCTS & INSTRUMENTATIONS

ARTHROVIMS® SOFT ANCHOR is packaged in a double pouch and sterilized by ethylene oxide. Single use. Do not re-sterilize. Instrumentation for arthroscopic or mini-open procedures is available. For any further information, please refer to the IFU.

DESIGNATION	SIZE	REFERENCE
ARTHROVIMS SOFT ANCHOR - 2 sutures	2,9	VIMS-029
ARTHROVIMS SOFT ANCHOR - 2 sutures with needles	2,9	VIMS-029N
ARTHROVIMS SOFT ANCHOR - 1 suture + 1 tape blue	2,9	VIMS-29FTB
ARTHROVIMS SOFT ANCHOR - 1 suture + 1 tape white/blue	2,9	VIMS-029FTWB
ARTHROVIMS SOFT ANCHOR - 2 sutures	1,8	VIMS-018
ARTHROVIMS SOFT ANCHOR - 2 sutures with needles	1,8	VIMS-018N
ARTHROVIMS SOFT ANCHOR - 1 suture with needle	1,5	VIMS-015N
ARTHROVIMS SOFT ANCHOR - 1 suture with needle	1,2	VIMS-012N

INSTRUMENTATION

ARTHROVIMS® SOFT ANCHOR must be used with the specific instrumentation provided by VIMS.

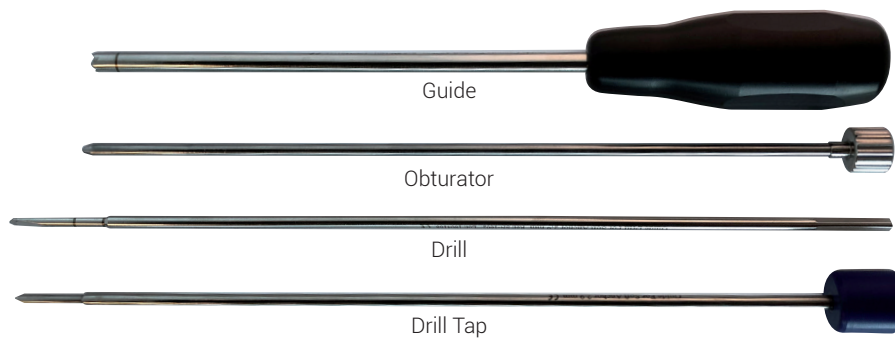

CE 2292
CLASSE IIb / ANCHOR

CE
CLASSE I / INSTRUMENTATION

Document design date : Février 2024 - Manufacturer : DORATEK - Range : All suture anchor - Brand Name : ARTHROVIMS SOFT ANCHOR, INSTRUMENTATION
User : Health professional - CE Mark N°: 2292 - Medical Device Class : IIb, I non stérile - Indications : Fixation of soft or prosthetic tissues to the bone - Recommendations for use: it is highly recommended to read the labels and the IFU.
BR_ARTHROVIMS_ANCHOR_TEXTILE_EN_V02



VIMS SA
10 AVENUE DE FONTRÉAL - EUROCENTRE
31620 VILLENEUVE LES BOULOC - FRANCE
TÉL. +33 5 34 45 09 09 - FAX. +33 5 61 22 85 84

VIMS-SYSTEM.COM

A ZIMMER BIOMET company